

Change 1/4 loop; N15 band filter, 15N-CIGAR

Sample Name:

Benzamide\_15N\_CIGAR

Data Collected on:

inovanmr.chem.wvu.edu-inova500

Archive directory:

/home/vnmr1/vnmrsys/data

Sample directory:

Benzamide\_15N\_CIGAR\_20140303\_01

FidFile: CIGAR\_01

Pulse Sequence: CIGAR

Solvent: dmsco

Data collected on: Mar 3 2014

Temp. 25.0 C / 298.1 K

Operator: vnmr1

Relax. delay 1.000 sec

Acq. time 0.150 sec

Width 8000.0 Hz

2D Width 2531.9 Hz

4 repetitions

2 x 200 increments

OBSERVE H1, 499.7522374 MHz

DATA PROCESSING


Sq. sine bell 0.075 sec

F1 DATA PROCESSING

Gauss apodization 0.006 sec

FT size 4096 x 2048

Total time 35 min

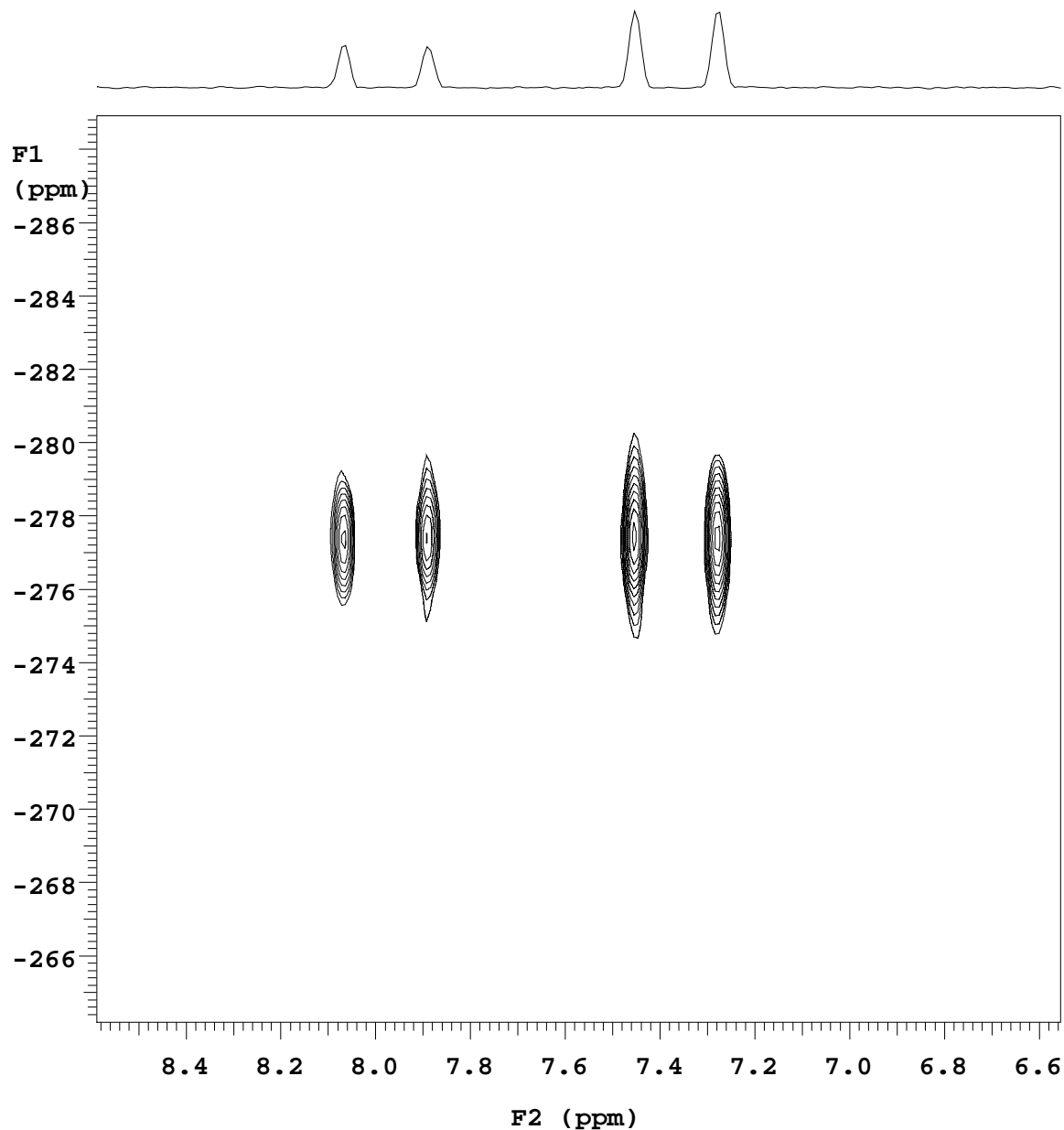

Plotname: CIGAR\_01\_plot01